




Activate NewCrop **WebView2** Interface

Following steps only once from any workstation

- 1) Go to PowerSoftMD website Support Tab, then Downloads button
- 2) Download and install Latest Updates to Main Computer

Following steps must be done at each workstation that will use eScripts


- 1) Go to PowerSoftMD website Support Tab, then Downloads button
- 2) Select Link: [Microsoft WebView2 Install](#)
- 3) Click and run the "Download" option under "Evergreen Bootstrapper"

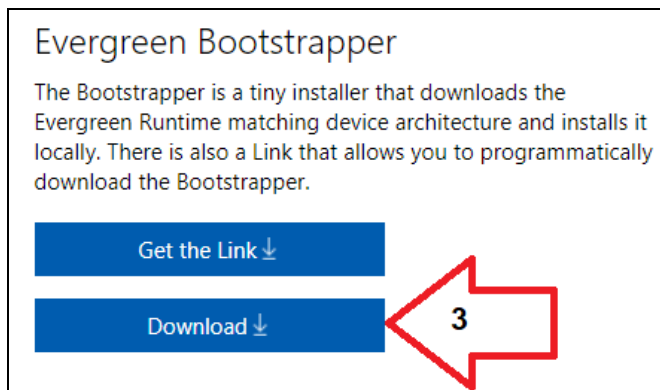


- 4) Select Link: [WebView2 Run Times Install](#)
Download and run to Drive C:
- 5) Run the PowerSoftMD workstation update **twice**

Following steps only once from any workstation


- 1) Login to PowerSoftMD, from the Primary Menu use the “Options” button

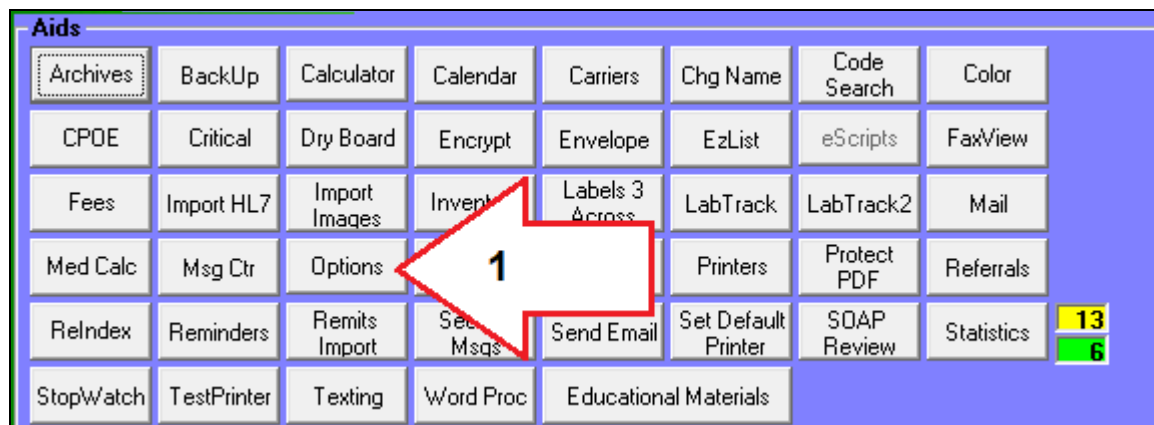

- 2) Check the “Use Browser CefSharp eScripts Interface”. Note: Even though you are setting up WebView2 Browser activate the CefSharp option! When prompted enter the Standard Special Password, then press “Save/Exit”

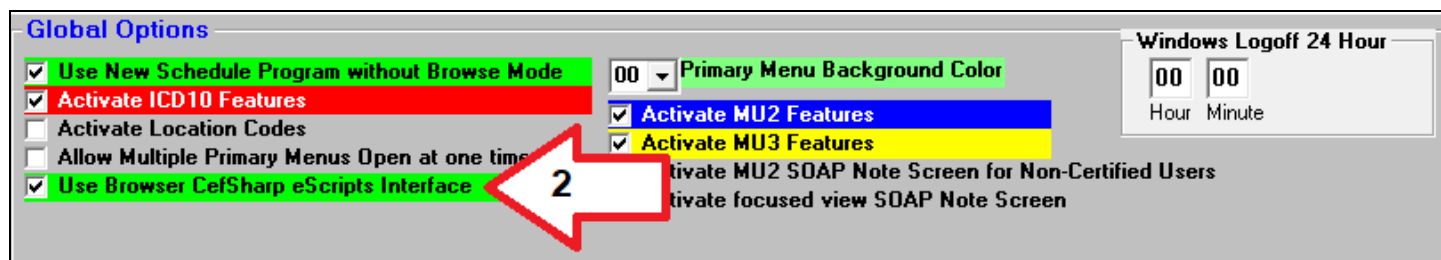


- 3) From any patient’s SOAP Note screen use the top toolbar “eScripts” option and then select “eScripts Menu”.

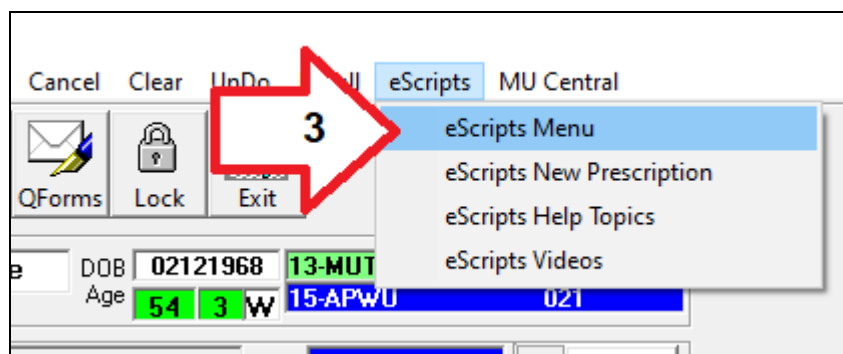

continued ...

4) on the “EdgeScripts” screen select “Tools”

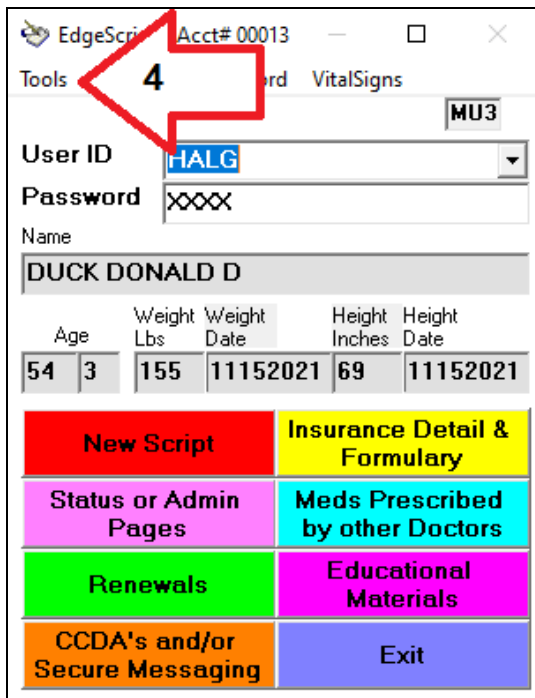

EdgeScripts: Acct# 00013

Tools **4** Help ChgPassword VitalSigns

MU3

User ID: HALG

Password: XXXX

Name: DUCK DONALD D

Age	Weight Lbs	Weight Date	Height Inches	Height Date	
54	3	155	11152021	69	11152021


New Script Insurance Detail & Formulary

Status or Admin Pages Meds Prescribed by other Doctors

Renewals Educational Materials

CCDA's and/or Secure Messaging Exit

5) On the drop down, list check the option “WebView2 Mode”, the password will be the the letter “E” plus the “Special Password”.

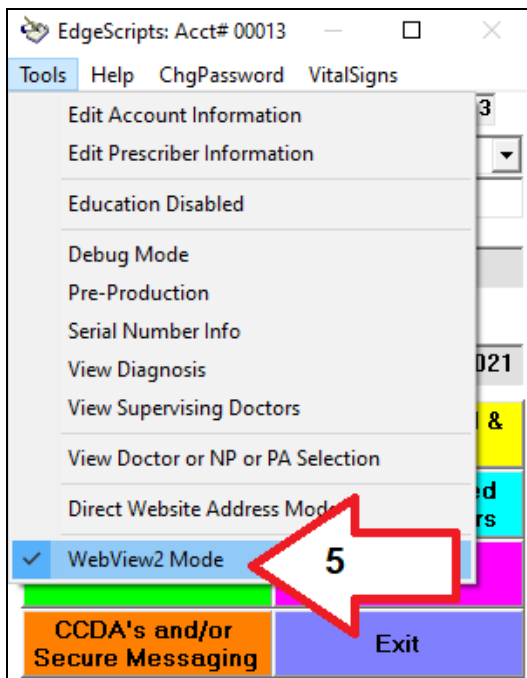

EdgeScripts: Acct# 00013

Tools Help ChgPassword VitalSigns

- Edit Account Information
- Edit Prescriber Information
- Education Disabled
- Debug Mode
- Pre-Production
- Serial Number Info
- View Diagnosis
- View Supervising Doctors
- View Doctor or NP or PA Selection
- Direct Website Address Mode
- ✓ WebView2 Mode 5**

CCDA's and/or Secure Messaging Exit